Directions for Use (Single Patient Use Only)

Product Reorder Numbers
Small VHC- PN 88810
Medium VHC- PN 78810

ENGLISH

Prior to use, carefully examine the product for proper assembly, foreign objects and damaged or missing parts. Ensure the mask adapter assembly is properly secured to the chamber body. Remove any foreign objects prior to use. The product should be replaced IMMEDIATELY if there are any damaged or missing parts. If necessary, use the Metered Dose Inhaler (MDI) alone until a replacement is obtained. If the patient's symptoms worsen, please seek immediate medical attention.

Remove cap from the MDI.
Before use, ensure the instructions supplied with the MDI have been read.

Shake the MDI immediately before each use as per the instructions supplied with MDI. Insert the MDI into the backpiece of the chamber.

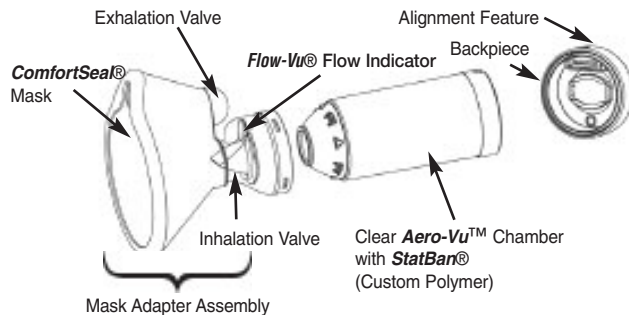
Apply mask to face and ensure there is an effective mask/face seal by watching the movement of the **Flow-Vu®** Inspiratory Flow Indicator. As the patient inhales through the chamber, the brightly colored flow indicator should move toward the patient, providing visual feedback that an effective mask/face seal has been obtained.

While maintaining an effective mask/face seal, the patient should exhale through the mask. The **Flow-Vu®** Inspiratory Flow Indicator will return to the vertical position when the patient stops inhalation or exhales. The exhaled breath exits the chamber through the Exhalation Valve mounted in the mask.

The MDI should be depressed as the patient begins an inhalation. Use the **Flow-Vu®** Inspiratory Flow Indicator to assist in the coordination of this maneuver. When the Flow Indicator starts to move toward the patient, this is the correct time to depress the MDI canister. Maintain the effective mask/face seal and ensure an additional 5-6 breaths are taken through the chamber before administering another puff of medication.

Ensure only one puff of medication is administered at a time and that the physician's dosage instructions are followed.

Patient should wait 30 seconds before administering another dose of medication (steps 3-5), or as instructed by their physician or health care provider.

Cleaning Instructions

ENGLISH

The **AeroChamber MAX®** VHC should be cleaned weekly, following the instructions provided below.

Remove the flexible backpiece and the mask adapter assembly as shown.

Do not remove the mask from the mask adapter assembly.

Soak the 3 parts for 15 minutes in a mild solution of liquid dish detergent and lukewarm clean water. Agitate gently.

Rinse parts in clean water.

Shake out excess water from the 3 parts and allow to air dry as shown.

Ensure parts are dry before reassembly.

To reassemble, fit the Mask Adapter assembly on the end of the chamber and twist firmly until components securely lock into position.

Center the alignment feature on the flexible backpiece with the flow indicator and exhalation valve as shown. Press firmly to seat backpiece.

Cautions

- Read all instructions before use.
- Do not leave **AeroChamber MAX®** VHC unattended with children.
- Do not disassemble the product beyond what is recommended in the cleaning instructions or damage may result.
- If unsure on how to use this product, talk to your physician or pharmacist or Health Care provider.
- The **AeroChamber MAX®** VHC is a medical device. To ensure proper performance it should only be cleaned according to these instructions. It should NOT be cleaned in a dishwasher.

Notes

- Storage and Operating Range: 40°F - 105°F at 15-95%R.H.
- VHC = Valved Holding Chamber
- Product may need to be replaced after 12 months of continuous use. Environmental conditions, storage and proper cleaning can affect product life span.

Indications for Use: The **AeroChamber MAX®** Valved Holding Chamber ("VHC") is intended to be used by patients who are under the care or treatment of a licensed health care professional or physician. The device is intended to be used by these patients to administer aerosolized medication from most pressurized Metered Dose Inhalers, prescribed by a physician or health care professional. The intended environments for use include the home, hospitals and clinics.

AeroChamber MAX®

Valved Holding Chamber with **Flow-Vu®** Indicator
and Small or Medium **ComfortSeal®** Mask

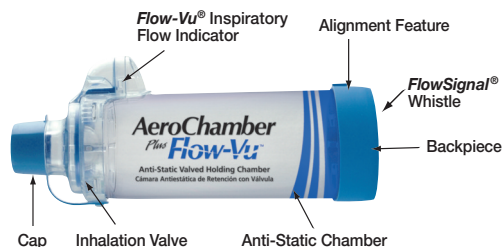
monaghan™

Manufactured by: Monaghan Medical Corporation
5 Latour Ave., Suite 1600 Plattsburgh, NY 12901
For additional information visit www.monaghanmed.com
Patent Numbers: 5,645,049; 5,848,588; 5,988,160; 6,293,279; 6,345,617; 6,435,177; 6,904,908; and patents pending
© 2004 Monaghan Medical Corporation
™ Trademarks and ® Registered Trademarks Printed in USA. P/N 88802-06 11/07

AeroChamber Plus Flow-Vu™

Anti-Static Valved Holding Chamber

MOUTHPIECE

MOUTHPIECE CHAMBER

INDICATIONS FOR USE

This product is intended to be used by patients who are under the care or treatment of a physician or licensed healthcare professional. The device is intended to be used by these patients to administer aerosolized medication from most pressurized Metered Dose Inhalers. The intended environments for use include the home, hospitals and clinics.

INSTRUCTIONS FOR USE

THIS PRODUCT CAN BE USED RIGHT OUT OF PACKAGE. BEFORE USE, ENSURE THESE INSTRUCTIONS AND THE INSTRUCTIONS SUPPLIED WITH THE METERED DOSE INHALER (MDI) HAVE BEEN READ.

<p>1</p> <p>Carefully examine the product for damage, missing parts or foreign objects. Remove any foreign objects prior to use. The product should be replaced IMMEDIATELY if there are any damaged or missing parts. If necessary, use the Metered Dose Inhaler (MDI) alone until a replacement is obtained. If the patient's symptoms worsen, please seek immediate medical attention.</p>	<p>2</p> <p>Remove cap(s) from the MDI and chamber (if applicable).</p>	<p>3</p> <p>Shake the MDI immediately before each use as per the instructions supplied with the MDI.</p>	<p>4</p> <p>Insert the MDI into the backpiece of the chamber.</p>	<p>5</p> <p>Put mouthpiece into mouth and close lips around it to ensure an effective seal. The Flow-Vu® Inspiratory Flow Indicator only moves if the patient has a good seal.</p>	<p>6</p> <p>Breathe out gently and depress the MDI at the beginning of a slow inhalation. Use the Flow-Vu® Indicator to assist in the coordination of this maneuver. Breathe in slowly and deeply through the mouth until a full breath has been taken. Hold breath for 5 – 10 seconds, if possible. Otherwise, keep lips tight on the mouthpiece to maintain seal for 2 – 3 breaths after the MDI is depressed. <i>Slow down inhalation if you hear the FlowSignal® Whistle sound.</i> Administer one (1) puff at a time.</p>	<p>7</p> <p>Follow instructions supplied with the MDI on how long to wait before repeating steps 3 – 6 as prescribed.</p>
---	--	---	--	--	--	--

CLEANING INSTRUCTIONS

THIS PRODUCT CAN BE USED RIGHT OUT OF THE PACKAGE AND THEN CLEANED WEEKLY.

<p>1</p> <p>Remove the backpiece only. Do not tamper with valve during cleaning or disassemble the product beyond what is recommended or damage may result.</p>	<p>2</p> <p>Soak the parts for 15 minutes in a mild solution of liquid dish detergent and lukewarm clean water. Agitate gently.</p>	<p>3</p> <p>Rinse parts in clean water.</p>	<p>4</p> <p>Shake out excess water from the parts and allow to air dry in a vertical position. Ensure parts are dry before reassembly.</p>	<p>5</p> <p>To reassemble, center the Alignment Feature on the backpiece with the Flow-Vu® Indicator, as shown. Press firmly to attach the backpiece. For mouthpiece models, the protective cap should always be placed on the mouthpiece when the product is not in use.</p>
--	--	--	---	---

Notes:

- Storage and operating range 5° C - 40° C (41° F - 104° F) at 15 to 95% relative humidity.
- Product may need to be replaced after 12 months of use. Environmental conditions, storage and proper cleaning can affect product life span.
- THIS PRODUCT CONTAINS NO LATEX.**
- Do not share this medical device.
- Clarity of the chamber is a result of the properties of the **StatBan®** anti-static material.
- If you notice medication build-up in your chamber, wash the inside of the chamber gently with a soft cloth.

Cautions:

- To ensure proper performance this product should only be cleaned according to these instructions.
- Do not leave the chamber unattended with children.
- Federal (USA) law restricts the sale of this device on or by the order of a physician.

Product Reorder Number:
79810/case 10, 79850/case 50

Rx Only

Mouthpiece

Manufactured by:

monaghan™

Monaghan Medical Corporation, 5 Latour Ave., Suite 1600, Plattsburgh, NY 12901
www.monaghanmed.com

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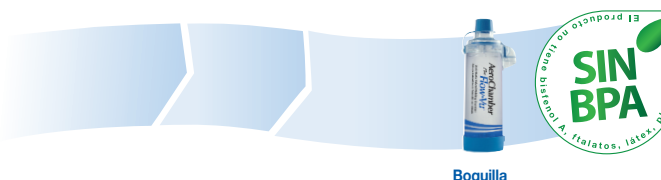
Boquilla

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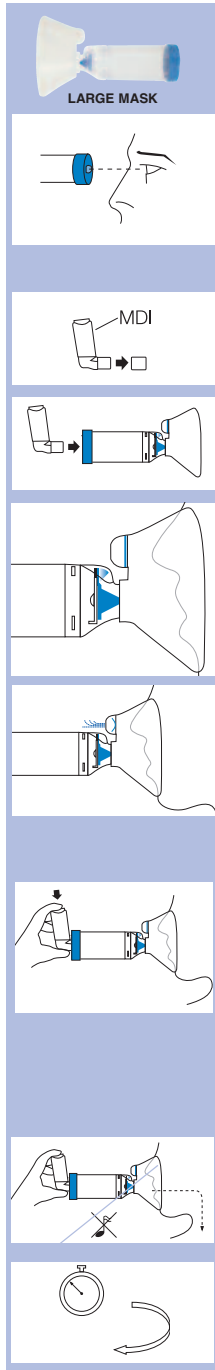
monaghan™

Monaghan Medical Corporation, 5 Latour Ave., Suite 1600, Plattsburgh, NY 12901
www.monaghanmed.com

Directions for Use (Single Patient Use Only)

Product Reorder Number
Large VHC- PN 80810

ENGLISH

Prior to use, carefully examine the product for proper assembly, foreign objects and damaged or missing parts. Ensure the mask adapter assembly is properly secured to the chamber body. Remove any foreign objects prior to use. The product should be replaced IMMEDIATELY if there are any damaged or missing parts. If necessary, use the Metered Dose Inhaler (MDI) alone until a replacement is obtained. If the patient's symptoms worsen, please seek immediate medical attention.

Remove cap from the MDI. Before use, ensure the instructions supplied with the MDI have been read.

Shake the MDI immediately before each use as per the instructions supplied with MDI. Insert the MDI into the back piece of the chamber.

Apply mask to face and ensure there is an effective mask/face seal by watching the movement of the **Flow-Vu®** Inspiratory Flow Indicator. As the patient inhales through the chamber, the brightly colored flow indicator should move toward the patient, providing visual feedback that an effective mask/face seal has been obtained.

While maintaining an effective mask/face seal, the patient should exhale through the mask. The **Flow-Vu®** Inspiratory Flow Indicator will return to the vertical position when the patient stops inhalation or exhales. The exhaled breath exits the chamber through the Exhalation Valve mounted in the mask.

The MDI should be depressed as the patient begins a slow inhalation. Use the **Flow-Vu®** Inspiratory Flow Indicator to assist in the coordination of this maneuver. When the Flow Indicator starts to move toward the patient, this is the correct time to depress the MDI canister. Maintain the effective mask/face seal and ensure an additional 2-3 breaths are taken through the chamber before administering another puff of medication.

Ensure only one puff of medication is administered at a time and that the physician's dosage instructions are followed.

Slow down inhalation if you hear the **FLOWSignal®** Whistle sound.

Follow instructions supplied with the MDI regarding the amount of time to wait before repeating instructions 3-6 as prescribed.

AeroChamber[®] MAX[®]

Valved Holding Chamber with **Flow-Vu®** Indicator
and Large **ComfortSeal®** Mask

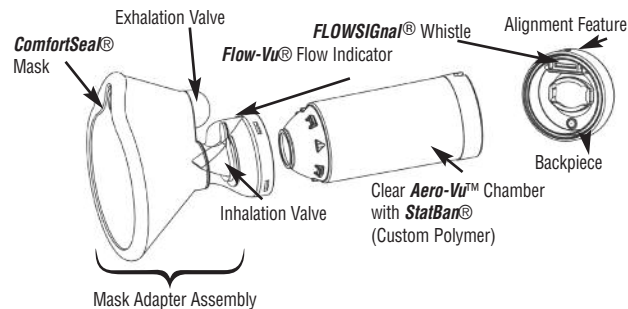
monaghan[™]

Manufactured by: Monaghan Medical Corp.
5 Latour Ave., Suite 1600 Plattsburgh, NY 12901
For additional information visit www.monaghanmed.com

© 2004 Monaghan Medical Corp.

™ Trademarks and ® Registered Trademarks. Printed in USA. P/N 80702-06 1/10


Cleaning Instructions

ENGLISH

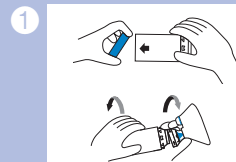The **AeroChamber MAX®** VHC should be cleaned weekly, following the instructions provided below.

Remove the flexible backpiece and the mask adapter assembly as shown.

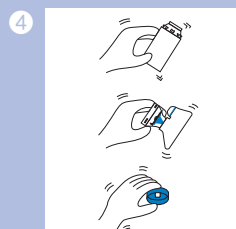
Do not remove the mask from the mask adapter assembly.


Soak the 3 parts for 15 minutes in a mild solution of liquid dish detergent and lukewarm clean water. Agitate gently.

Rinse parts in clean water.

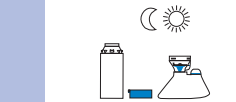
Shake out excess water from the 3 parts and allow to air dry as shown.

Ensure parts are dry before reassembly.

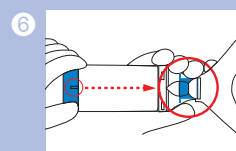To reassemble, fit the Mask Adapter assembly on the end of the chamber and twist firmly until components securely lock into position.

Center the alignment feature on the flexible backpiece with the flow indicator and exhalation valve as shown. Press firmly to seat backpiece.

Cautions

- Read all instructions before use.
- Do not leave **AeroChamber MAX®** VHC unattended with children.
- Do not disassemble the product beyond what is recommended in the cleaning instructions or damage may result.
- If unsure on how to use this product, talk to your physician or pharmacist or Health Care provider.
- The **AeroChamber MAX®** VHC is a medical device. To ensure proper performance it should only be cleaned according to these instructions. It should NOT be cleaned in a dishwasher.

Notes

- Storage and Operating Range: 40°F - 105°F at 15-95%R.H.
- VHC = Valved Holding Chamber
- Product may need to be replaced after 12 months of continuous use. Environmental conditions, storage and proper cleaning can affect product life span.

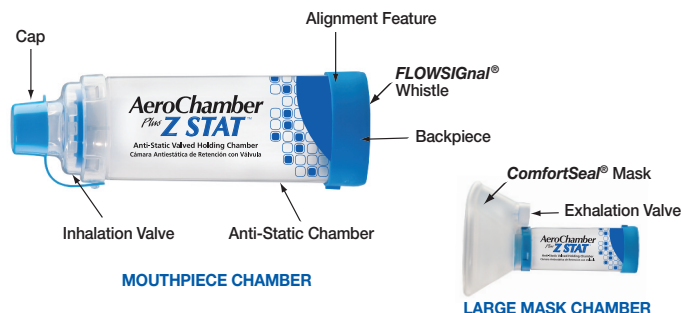
Indications for Use: The **AeroChamber MAX®** Valved Holding Chamber ("VHC") is intended to be used by patients who are under the care or treatment of a licensed health care professional or physician. The device is intended to be used by these patients to administer aerosolized medication from most pressurized Metered Dose Inhalers, prescribed by a physician or health care professional. The intended environments for use include the home, hospitals and clinics.

Rx Only
LATEX FREE

AeroChamber Plus Z STAT™

Anti-Static Valved Holding Chamber

MOUTHPIECE
LARGE MASK

INDICATIONS FOR USE

This product is intended to be used by patients who are under the care or treatment of a physician or licensed healthcare professional. The device is intended to be used by these patients to administer aerosolized medication from most pressurized Metered Dose Inhalers. The intended environments for use include the home, hospitals and clinics.

INSTRUCTIONS FOR USE

THIS PRODUCT CAN BE USED RIGHT OUT OF PACKAGE. BEFORE USE, ENSURE THESE INSTRUCTIONS AND THE INSTRUCTIONS SUPPLIED WITH THE METERED DOSE INHALER (MDI) HAVE BEEN READ.

MOUTHPIECE	1	2	3	4	5	6	7
	<p>Carefully examine the product for damage, missing parts or foreign objects. Remove any foreign objects prior to use. The product should be replaced IMMEDIATELY if there are any damaged or missing parts. If necessary, use the Metered Dose Inhaler (MDI) alone until a replacement is obtained. If the patient's symptoms worsen, please seek immediate medical attention.</p>	<p>Remove cap(s) from the MDI and chamber (if applicable).</p>	<p>Shake the MDI immediately before each use as per the instructions supplied with the MDI.</p>	<p>Insert the MDI into the backpiece of the chamber.</p>	<p>Put mouthpiece into mouth and close lips around it to ensure an effective seal.</p>	<p>Breathe out gently and depress the MDI at the beginning of a slow inhalation. Breathe in slowly and deeply through the mouth until a full breath has been taken. Hold breath for 5 – 10 seconds, if possible. Otherwise, keep lips tight on the mouthpiece to maintain seal for 2 – 3 breaths after the MDI is depressed. <i>Slow down inhalation if you hear the FLOWSIGNAL® Whistle sound. Administer one (1) puff at a time.</i></p>	<p>Follow instructions supplied with the MDI on how long to wait before repeating steps 3 – 6 as prescribed.</p>
LARGE MASK	1	2	3	4	5	6	7
	<p>Carefully examine the product for damage, missing parts or foreign objects. Remove any foreign objects prior to use. The product should be replaced IMMEDIATELY if there are any damaged or missing parts. If necessary, use the Metered Dose Inhaler (MDI) alone until a replacement is obtained. If the patient's symptoms worsen, please seek immediate medical attention.</p>	<p>Remove cap(s) from the MDI and chamber (if applicable).</p>	<p>Shake the MDI immediately before each use as per the instructions supplied with the MDI.</p>	<p>Insert the MDI into the backpiece of the chamber.</p>	<p>Apply mask to face and ensure an effective seal.</p>	<p>Breathe out gently and depress the MDI at the beginning of a slow inhalation. Maintain seal for 5 – 6 breaths after the MDI is depressed. <i>Slow down inhalation if you hear the FLOWSIGNAL® Whistle sound. Administer one (1) puff at a time.</i></p>	<p>Follow instructions supplied with the MDI on how long to wait before repeating steps 3 – 6 as prescribed.</p>

CLEANING INSTRUCTIONS FOR MASK AND MOUTHPIECE CHAMBERS

THIS PRODUCT CAN BE USED RIGHT OUT OF THE PACKAGE AND THEN CLEANED WEEKLY.

1	2	3	4	5
<p>Remove the backpiece only. Do not tamper with valve during cleaning or disassemble the product beyond what is recommended or damage may result. For mask product, do not remove mask.</p>	<p>Soak the parts for 15 minutes in a mild solution of liquid dish detergent and lukewarm clean water. Agitate gently.</p>	<p>Rinse parts in clean water.</p>	<p>Shake out excess water from the parts and allow to air dry in a vertical position. Ensure parts are dry before reassembly.</p>	<p>To reassemble, press firmly to attach the backpiece. For mouthpiece models, the protective cap should always be placed on the mouthpiece when the product is not in use.</p>

Notes:

- Storage and operating range 5° C - 40° C (41° F - 104° F) at 15 to 95% relative humidity.
- Product may need to be replaced after 12 months of use. Environmental conditions, storage and proper cleaning can affect product life span.
- THIS PRODUCT CONTAINS NO LATEX.
- Do not share this medical device.
- Clarity of the chamber is a result of the properties of the StatBan® anti-static material.
- If you notice medication build-up in your chamber, wash the inside of the chamber gently with a soft cloth.

Cautions:

- To ensure proper performance this product should only be cleaned according to these instructions.
- Do not leave the chamber unattended with children.
- Federal (USA) law restricts the sale of this device on or by the order of a physician.

Product Reorder Number:

MOUTHPIECE: 79710/case 10, 79750/case 50

LARGE MASK: 80710/case 10

Rx Only

© 2005, 2010 Monaghan Medical Corporation™ and ® are trademarks and registered trademarks of Monaghan Medical Corporation or an affiliate of Monaghan Medical Corporation. This device was Made in the USA from US and Canadian parts.

Manufactured by:

monaghan™

Monaghan Medical Corporation, 5 Latour Ave., Suite 1600, Plattsburgh, NY 12901
www.monaghanmed.com

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ENGLISH

Carefully examine the product for damage, missing parts or foreign objects. Any foreign objects should be removed. The product should be replaced **IMMEDIATELY** if there are any damaged or missing parts. If necessary, use the inhaler (MDI) alone until a replacement is obtained. If your symptoms worsen, please seek immediate medical attention.

Before use, make sure that instructions supplied with inhaler (MDI) have been read. Remove cap(s).

Shake the inhaler (MDI) well immediately before each use. Insert the inhaler (MDI) into the back piece of the chamber.

Put mouthpiece into mouth.

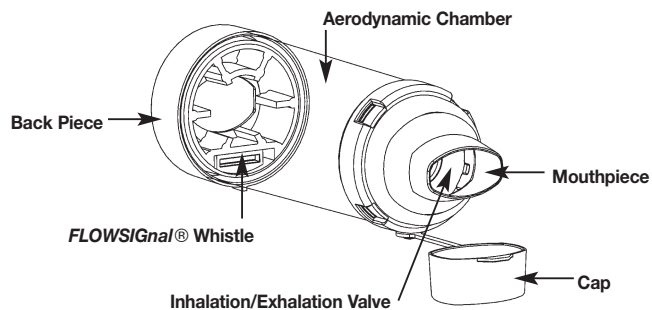
Breathe out gently and depress inhaler (MDI) at beginning of slow deep inhalation. Hold breath as long as possible, up to 10 seconds before breathing out.

If patient has difficulty with slow deep breaths, an alternative is to keep mouth tight on mouthpiece and breathe slowly 2-3 times after depressing inhaler (MDI).

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Follow instructions supplied with the inhaler (MDI) on amount of time to wait before repeating steps 3–6 as prescribed.

Cleaning Instructions

ENGLISH

Clean **AeroChamber Plus®** VHC only as per instructions before first use, then weekly.

Remove back piece only. Do not tamper with valve during cleaning.

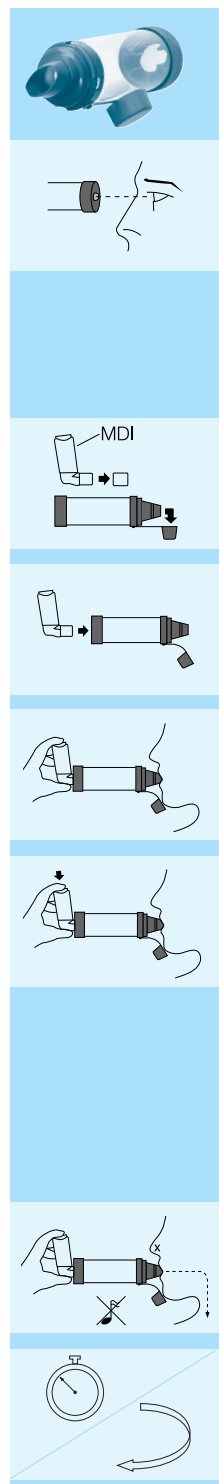
Soak both parts for 15 minutes in lukewarm water with liquid detergent. Agitate gently.

Rinse mouthpiece only in water.

Shake out excess water.
Do not rub dry.

Let air dry in vertical position.

Replace back piece when unit is completely dry and ready for use.

AeroChamber® Plus®

Valved Holding Chamber with **FLOWSignal®** Whistle

monaghan™

Monaghan Medical Corporation, 5 Latour Ave., Suite 1600, Plattsburgh, NY 12901
www.monaghanmed.com
 © 1999, 2002 Monaghan Medical Corporation
 ™ trademarks and ® registered trademarks
 Covered by one or more of the following Patents
 5,848,588; 6,345,617; 6,293,279; 6,435,177; and patents pending
 Part Number: 76071 AWC34 Revision: 10/09

Cautions

- Read all instructions before use.
- Do not leave **AeroChamber Plus®** VHC unattended with children. This is a medical device, not a toy.
- Do not disassemble the product beyond what is recommended in the Cleaning Instructions or damage may result.
- If unsure on how to use this product, talk to your physician or pharmacist.

Notes

- VHC = Valved Holding Chamber
- **This product contains no latex.**
- The **AeroChamber Plus®** VHC may need to be replaced after 12 months of continuous use. Environmental conditions, storage and proper cleaning may affect life span.

Indications for Use: The **AeroChamber Plus®** Valved Holding Chamber ("VHC") is intended to be used by patients who are under the care or treatment of a licensed health care professional or physician. The device is intended to be used by these patients to administer aerosolized medication from most pressurized Metered Dose Inhalers, prescribed by a physician or health care professional. The intended environments for use include the home, hospitals and clinics.

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<p>1</p>	<p>2</p>	<p>3</p>	<p>4</p>	<p>5</p>	<p>6</p>	<p>7</p>
<p>Carefully examine the product for damage, missing parts or foreign objects. Remove any foreign objects prior to use. The product should be replaced IMMEDIATELY if there are any damaged or missing parts. If necessary, use the Metered Dose Inhaler (MDI) alone until a replacement is obtained. If the patient's symptoms worsen, please seek immediate medical attention.</p>	<p>Remove cap from the MDI.</p>	<p>Shake the MDI immediately before each use as per the instructions supplied with the MDI.</p>	<p>Insert the MDI into the backpiece of the chamber.</p>	<p>Apply mask to face and ensure an effective seal.</p>	<p>Breathe out gently and depress the MDI at the beginning of a slow inhalation. Maintain seal for 5 - 6 breaths after the MDI is depressed. Administer one (1) puff at a time.</p>	<p>Follow instructions supplied with the MDI on how long to wait before repeating steps 3 - 6 as prescribed.</p>

TIPS

- 1 Give praise and rewards.
- 2 Use teddy bears or dolls to 'pretend' to use the chamber to make it look like a game.
- 3 If possible, give treatments when the child is happy and not crying.

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<p>1</p>	<p>2</p>	<p>3</p>	<p>4</p>	<p>5</p>
<p>Remove the backpiece only. Do not tamper with valve during cleaning or disassemble the product beyond what is recommended or damage may result.</p>	<p>Soak the parts for 15 minutes in a mild solution of liquid dish detergent and lukewarm clean water. Agitate gently.</p>	<p>Rinse parts in clean water.</p>	<p>Shake out excess water from the parts and allow to air dry in a vertical position. Ensure parts are dry before reassembly.</p>	<p>To reassemble, press firmly to attach the backpiece.</p>

Notes:

- Storage and operating range 5° C - 40° C (41° F - 104° F) at 15 to 95% relative humidity.
- Product may need to be replaced after 12 months of use. Environmental conditions, storage and proper cleaning can affect product life span.
- THIS PRODUCT CONTAINS NO LATEX.**
- Do not share this medical device.
- Clarity of the chamber is a result of the properties of the StatBan® anti-static material.
- If you notice medication build-up in your chamber, wash the inside of the chamber gently with a soft cloth.

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Product Reorder Number:
SMALL MASK: 88710/case 10
MEDIUM MASK: 78710/case 10

Rx Only

© 2005, 2010 Monaghan Medical Corporation™ and ® are trademarks and registered trademarks of Monaghan Medical Corporation or an affiliate of Monaghan Medical Corporation. This device was Made in the USA from Canadian and US parts.

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