

DIRECTIONS FOR USE



FIRST DOSE ADMINISTRATION

Pull to remove **GREEN** caps labeled "1" and "2."

Place rounded tip against the middle of the outer thigh. The needle may be injected through clothing. Press down hard until auto-injector fires. Hold in place while slowly counting to 10.


Immediately call 911.

Never place thumbs, fingers, or hands over rounded tip.

Refer to Patient Information included with your Twinject device for more information. If you have questions, please contact 1-888-TWINJCT.

SECOND DOSE ADMINISTRATION

If symptoms have not improved within about 10 minutes after the first injection, a second dose is needed.

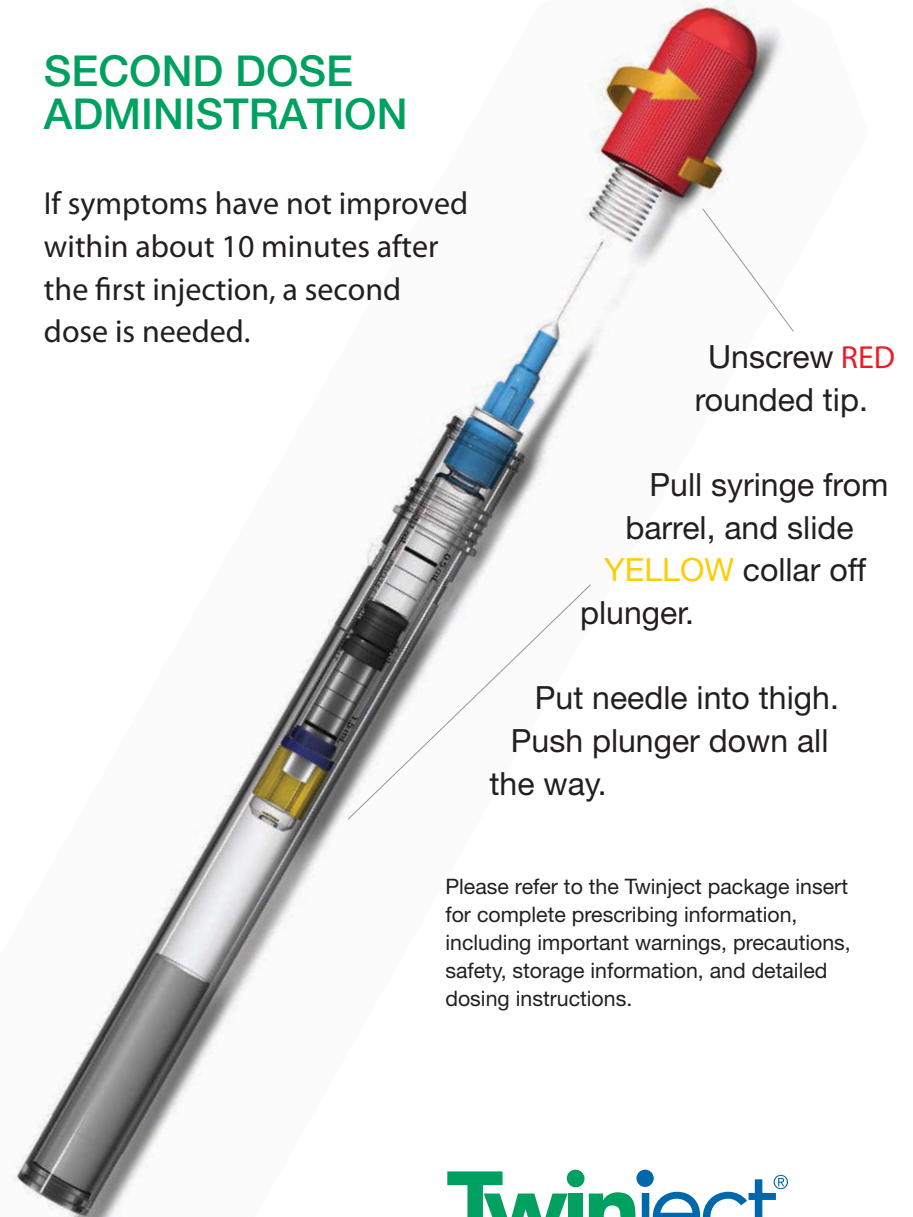

Unscrew **RED** rounded tip.

Pull syringe from barrel, and slide **YELLOW** collar off plunger.

Put needle into thigh. Push plunger down all the way.

Please refer to the Twinject package insert for complete prescribing information, including important warnings, precautions, safety, storage information, and detailed dosing instructions.

Twinject[®]
auto-injector
(epinephrine injection, USP 1:1000)

Available as: **0.3 mg** **0.15 mg**

Twinject[®] is a registered trademark of Sciele Pharma.
© 2008 Sciele Pharma, Inc., Atlanta, Georgia
All rights reserved.